

LABORATORY WORK SHEET

Name of the Student : Mohammad Mateen

Class : CSE - C Semester : 01

Course Code : A11SD03 Course Name : Chemistry lab

Name of the Course Faculty : Dr VNSR Venkateshwarrao

Faculty ID : JARE 10682

Exercise Number : 06

Week Number : 06

Date : 14/12/23

Roll Number									
2	3	9	5	1	A	0	5	5	B

DAY TO DAY EVALUATION:

Marks	Aim / Preparation	Algorithm / Procedure	Source Code	Program Execution	Viva - Voce	Total
		Performance in the Lab	Calculations and Graphs	Results and Error Analysis		
Max. Marks	4	4	4	4	4	20
Obtained	4	4	4	4	4	20

Signature of Faculty

START WRITING FROM HERE :


Estimation of Fe^{+2} by potentiometry Using $K_2Cr_2O_7$

Aim: To estimate the amount of ferrous iron present in the whole of the given solution by potentiometry.

Principle: * This is an example of redox titration and is based on the oxidation-reduction reaction between the titrand and the titrant. Here the end point is detected using a potentiometer.


* Potassium dichromate is an oxidising agent and in acid medium; it follows the half reaction to give $Cr(III)$ as the reduction product.

The total reaction is

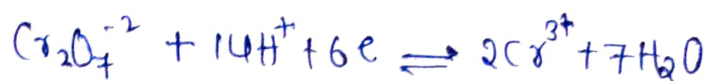

while Fe^{2+} which is used to titrate $\text{K}_2\text{Cr}_2\text{O}_7$ gets oxidised to Fe^{3+} as per the reaction. $\text{Fe}^{2+} \rightleftharpoons \text{Fe}^{3+} + e^-$

The overall ionic equation of this titration can be obtained by adding the above two:

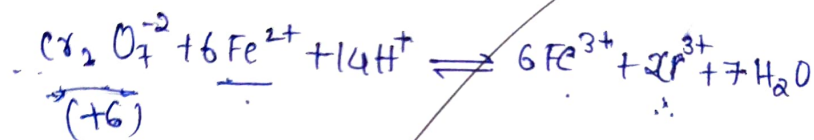

Preparation of Solutions

ii) Preparation of standard $K_2Cr_2O_7$ and ammonium iron(II) Sulphate

Using 100cm³ volumetric flasks prepare of 0.02M potassium dichromate solution and 0.10M ammonium iron(II) sulphate solution. You may have to add sufficient amount of dilute acid to prepare ammonium iron (II) sulphate solution.

2) Preparation of $2N H_2SO_4$ solution:

56ml of conc. H_2SO_4 is added Drop by Drop to 1000ml of water in a beaker by keeping it in a trough of water, this reaction is highly exothermic.

Note : H_2SO_4 must not be added to water.

Procedure:

Potentiometric Titration: *The given Fe^{2+} solution is made up to the mark of the given 100ml standard flask and shake the flask well for uniform concentration.

* pipette out 20ml of the Fe^{+2} solution into a clean 250ml beaker and add equal volume (20ml) of dil H_2SO_4 along with 100ml of distilled water to enable the electrodes to immerse well in the solution.


* A platinum electrode (indicator electrode) and a standard calomel electrode (reference electrode) from the potentiometer are dipped into the beaker. The solution in the beaker is stirred using a magnetic stirrer. The initial EMF is noted. 0.5ml of $\text{K}_2\text{Cr}_2\text{O}_7$ solution is added from the burette at regular intervals of time, while stirring the solution and the EMF is measured. The volume of $\text{K}_2\text{Cr}_2\text{O}_7$ added and the corresponding EMF readings are noted. At the end point there is a sharp increase in EMF due to the complete oxidation of the Fe^{+2} to Fe^{+3} . The addition of $\text{K}_2\text{Cr}_2\text{O}_7$ is continued till the equivalent point is crossed by at least 5ml.

S.NO.	Volume of $\text{K}_2\text{Cr}_2\text{O}_7$ added	EMF
1	0	0.301
2	1	0.360
3	2	0.377
4	3	0.392
5	4	0.403
6	5	0.413
7	6	0.423
8	7	0.433
9	8	0.446
10	9	0.463
11	10	0.496

12
13
14
15
16
17
18
19
20
21

11
12
13
14
15
16
17
18
19
20

0.641 → ↑↑
0.697
0.698
0.711
0.721
0.729
0.787
0.740
0.743
0.745

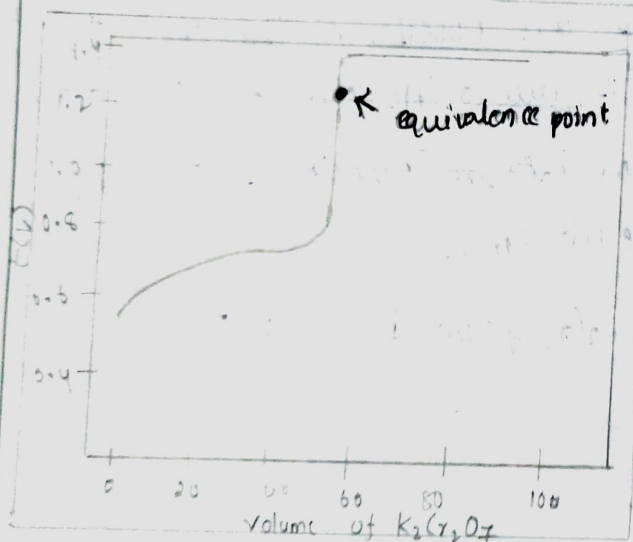

Modal graph

calculations:

$$V_1 N_1 = V_2 N_2$$

$$V_1 = \text{Volume of } K_2Cr_2O_7 = 10.5 \text{ ml}$$

$$N_1 = \text{Normality of } K_2Cr_2O_7 = 0.1 \text{ N}$$

$$V_2 = \text{Volume of } Fe^{+2} = 10 \text{ ml}$$

$$N_2 = \text{Normality of } Fe^{+2} = ?$$

$$N_2 = 4 \times (20 \times \text{Normality of } K_2Cr_2O_7) / 20$$

$$N_1 V_1 = N_2 V_2$$

$$0.1 \times 10.5 = N_2 \times 10$$


$$N_2 = \frac{10.5 \times 0.1}{10} = \frac{1.05}{10} = 0.105 \text{ N}$$

Result: Amount of Fe^{+2} present in the given solution = $N_2 \times \text{equivalent wt}$
 $= N_2 \times 55.85 \frac{\text{gms}}{100 \text{ ml}}$

⑥ Estimation of Fe^{+2} By Potentiometry Using $\text{K}_2\text{Cr}_2\text{O}_7$ scale

on x axis \rightarrow (units = ml)
on y axis \rightarrow 1 unit = 0.01V

\updownarrow
(V) $\times 100$

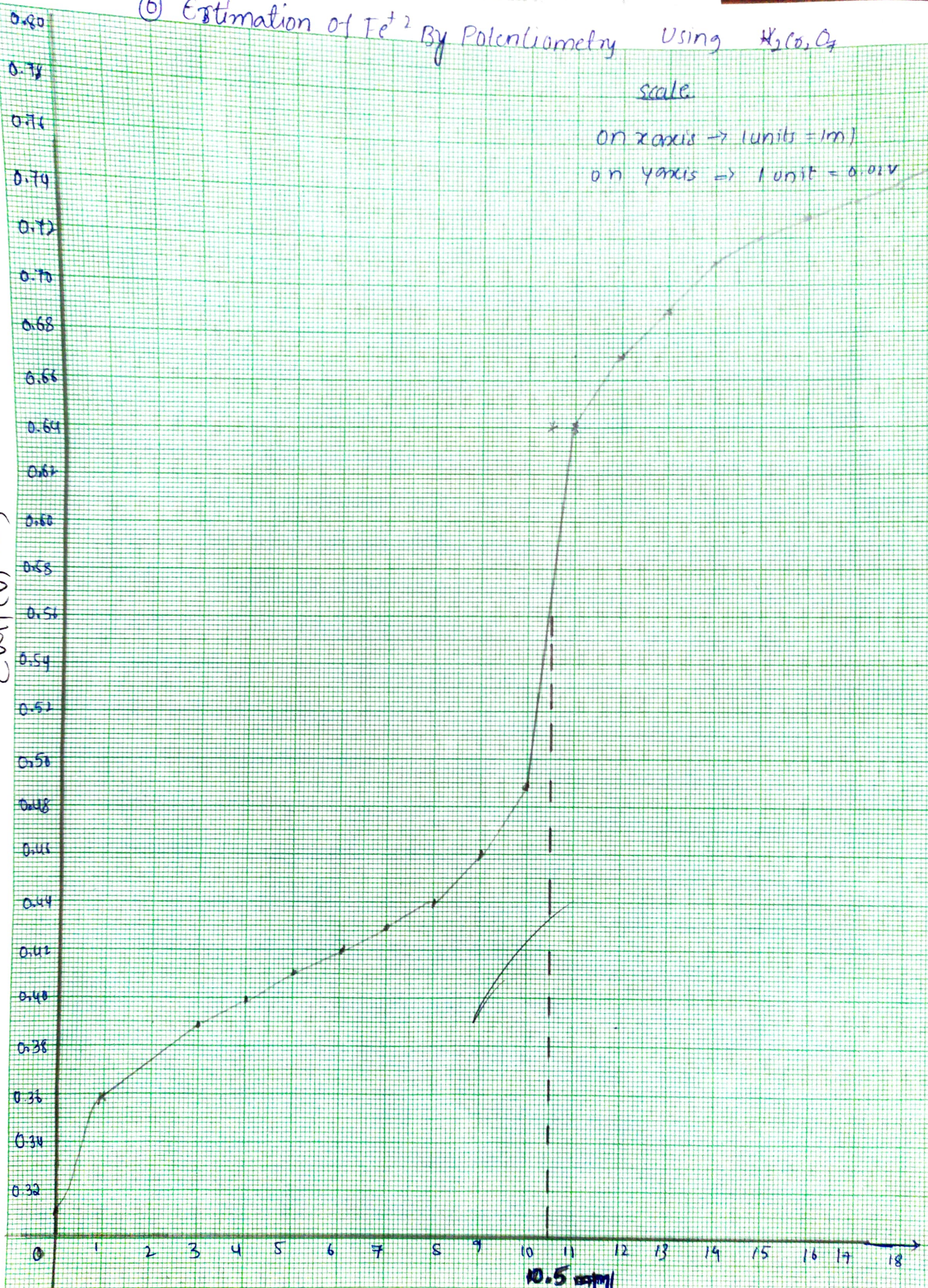

Volume of $\text{K}_2\text{Cr}_2\text{O}_7$ (ml) \rightarrow